# AquaFloss Station O300 Ozone

---

**EN - USER MANUAL**

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**HR - UPUTE ZA UPOTREBU**

---


## EN - USER MANUAL

---

### BOX CONTENTS

- TrueLife AquaFloss Station O300 Ozone
- 3x Standard nozzles, 1x Dental Plaque nozzle, 1x Orthodontic nozzle, 1x Periodontal nozzle
- Battery cover opener
- 1x AAA battery
- Manual and safety instructions

### PRODUCT DESCRIPTION

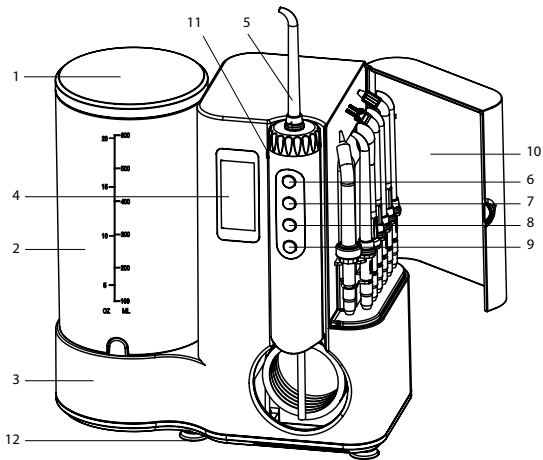
1. Water reservoir cover
2. Water reservoir
3. Base
4. Display
5. Nozzle
6. Power button

7. Increase power button
8. Decrease power button
9. Oxidation on/off button
10. Storage compartment door
11. Nozzle release button
12. Rubber suction cups

## Nozzles

**Please note:** There may be residual water in the device. This is distilled water that is used when testing the device.
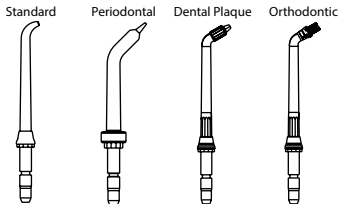
## INSTRUCTIONS FOR USE

Remove all packaging materials before first use.

If the water flosser has been exposed to a low temperature environment during transport, leave it at room temperature for at least two hours before use.

1. Connect the nozzle to the handle. It will click when correctly inserted.
2. Remove the water reservoir from the base and fill it with water not exceeding 40°C. The capacity of the reservoir is 600 ml. Return the reservoir to the base and make sure that no water overflows.
3. Plug in the water flosser.
4. Use the power increase and decrease buttons to set the pressure you want. The currently set pressure is shown on the display with a choice of 10 levels.
5. Use the oxidation on/off button to choose which water you want to use for cleaning. The „O3 Water“ icon will flash on the display for oxidised water and the „Natural Water“ icon for conventional water.
6. Lean over the sink and put the nozzle in your mouth. Aim the nozzle at your teeth and use the power button to start cleaning.
7. Leave your mouth slightly open to allow water to flow freely from your mouth into the sink. Direct the water stream at a 90° angle towards the teeth and gums. Slowly move the water stream around the teeth. Direct the water stream towards the gums and the spaces between the teeth. If you have braces, you can direct the water between the joints of the braces.
8. The water flosser has a preset automatic cleaning timer of 3 minutes. If you want to floss your teeth for a longer time, restart with the power button.
9. When you have finished cleaning, switch the device off using the power button. Open the lid of the water reservoir and pour out the remaining water. Switch the device on again to empty the water from the pump.

**Note:** To prevent the spread of bacteria, make sure that no water is left in the device. If you are not going to use the device for a long period of time, make sure it is thoroughly dried before storing it.

## INSERTING/REPLACING THE NOZZLE

Insert the nozzle into the opening at the top of the flosser. When correctly inserted you'll hear a click.

To remove the nozzle from the handle, press and hold the nozzle release button and pull the nozzle out.

The nozzles can be swivelled 360°.

The box includes the following nozzles:

- **Standard Jet** – used to rinse hard-to-reach areas and places where a conventional toothbrush cannot reach.
- **Orthodontic Jet** – a special orthodontic nozzle designed to clean and rinse hard-to-reach areas around braces.
- **Periodontal Jet** – directs the stream of water or mouthwash below the gum line thanks to the tapered silicone end.
- **Dental Plaque Jet** – used for more effective removal of dental plaque with easier access to affected areas.

**Note:** Only insert and remove the nozzles when the flosser is switched off.

**Note:** Nozzles are subject to personal hygiene and are not subject to warranty.

## REPLACING THE BATTERY IN THE HANDLE

The handle is powered by an AAA battery. When the controls stop working, it means the battery is dead and needs to be replaced.

Follow the instructions below:

1. Use the special opener to open the cover on the bottom of the handle.
2. Replace the battery with a new one.
3. Close and secure using the opener again.

We recommend removing the battery from the handle if you are not using the device for an extended period of time.

## MAINTENANCE

Use only water and neutral detergents for maintenance purposes. Never use caustic or abrasive cleaners as they can damage the equipment. Do not use water at temperatures above 50°C. Keep the equipment away from high temperatures and avoid direct sunlight.

### 1. Base

Remove dirt by wiping it down with a towel or cloth with liquid soap and water. Then wipe it dry. Do not immerse the device in water.

### 2. Water reservoir

Wash the reservoir with water and wipe with a dry cloth.

Dry the reservoir thoroughly if you are not going to use it for a long time.

After using a special product in the reservoir (mouthwash), rinse the reservoir each time. Partially fill the reservoir with clean warm water, direct the nozzle into the sink, turn on and let it flow until the reservoir is completely empty. Regular flushing prevents clogging and prolongs the life of the product.

### 3. Nozzles

Wash the nozzle with water and wipe with a cloth.

## TROUBLESHOOTING

<b>Problem</b>	<b>Cause</b>	<b>Solution</b>
Water leaks at the connection between the nozzle and the handle	The nozzle was incorrectly inserted	Remove nozzle and reattach
Periodontic nozzle is damaged	Nozzle is old	Replace the nozzle with new one
Water leaks from the hose/handle	Hose/handle is damaged	Return the device
The reservoir leaks	Valve on the reservoir is missing or incorrectly inserted	Install the valve correctly. If there is no valve, contact your dealer
Low water pressure	The reservoir is not seated correctly	Make sure that the outlet is working. Check the functionality of the electrical socket with other appliances
Water is leaking from the device	The motor is switched off	Empty the reservoir and the residual water from the hose.
The water smells	Dirt has settled in the nozzle, reservoir or hose	Wash the reservoir, rinse the hose with clean warm water and replace the nozzle

## SPECIFICATIONS

<b>Power consumption</b>	18 W
<b>Voltage</b>	AC100~240 V
<b>Frequency</b>	50/60 Hz
<b>Water resistance</b>	IPx7
<b>Pressure range</b>	30-125 PSI
<b>Reservoir volume</b>	600 ml
<b>Weight</b>	955 g
<b>Dimensions</b>	192×115×271 mm

## PRODUCT LIABILITY

TrueLife products are guaranteed to be free from manufacturing defects. In the event of such defects, elem6 s.r.o. will repair or replace the defective part or product with a comparable part or product at its sole discretion. The warranty will be deemed void if the device is altered or tampered with. The warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss, or other claims related to or resulting from the use of this product. elem6 s.r.o. shall in no event be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Return or exchange of spare parts and / or products may be subject to shipping, handling, and replacement or replenishment fees.

This product should be stored out of the reach of children. Use by children should take place only under the direct supervision of an adult.

## SAFETY GUIDELINES WHEN OPERATING THE DEVICE

1. Read the enclosed safety instructions carefully.
2. The latest version of the manual, including instructions for operating the device, can be downloaded at: [www.truelife.eu](http://www.truelife.eu).
3. Printing errors and changes of information reserved.

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