# truelife


## FitScale W4 BT

### EN - INSTRUCTIONS FOR USE

The TrueLife FitScale W4 BT wireless smart scale measures body weight using bioelectrical impedance analysis (BIA). BIA technology is used to measure body fat and obtain data on the total body composition. The TrueLife FitScale W4 BT measures not only weight and body fat ratio, but also fat free mass, subcutaneous fat, visceral fat, muscle mass, skeletal muscle, bone mass, water, protein, BMI, AMR and BMR calorie analysis, metabolic age and body type. The scale is intended for home use only.

**GETTING STARTED:**  
– Before turning on the device, insert the two included AAA 1.5V batteries into the scale according to the diagram under the battery cover.  
– Place the scale on a flat and stable surface or mat.

**Note: If you use the scale on a soft surface or carpet, it will not show the results correctly.**

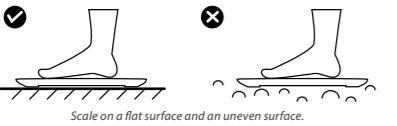Scale on a flat surface and an uneven surface.


– To achieve consistent and accurate measurement results, use the scale at the same time of day under approximately the same conditions. If possible, avoid weighing yourself in the following cases:

1. During or immediately after strenuous exercise;
2. Immediately after a shower, bath or sauna;
3. After heavy alcohol consumption;
4. During or immediately after a meal;
5. During a fever.

### INSTALLING THE APP AND WEIGHING YOURSELF FOR THE FIRST TIME:

– The TrueLife FitScale W4 BT scale communicates with your smartphone using Bluetooth wireless technology. For the proper operation of the application, you will need a smartphone running iOS 8.0 or higher, or Android 4.4 or higher.

#### 1. Download the app and register or log in:

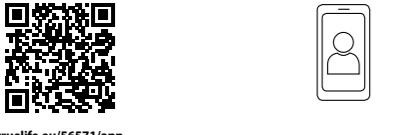1. Scan the QR code or search for the app in App Store or on Google Play  
2. Register and fill in your personal data and parameters.

#### 2. Pair the scale with the application:

1. First, turn on location and Bluetooth on your smartphone.  
2. Stand on the scales. The scale switches on automatically when stepping on to it.  
3. Then find the scale in the application and pair it.

#### 3. Weigh yourself:

- In order for the scale to remember you and pair it with your account, you need to weigh yourself immediately after registration. When weighing yourself for the first time, CODE 1 appears on the top left of the screen.
- Weight units can be switched from kilograms (kg) to pounds (lb) by pressing the button on the bottom of the scale. The units can also be changed in the application.


**Note: Each time you weigh yourself, first touch the scale so that the display shows 0.00 kg. Only then weigh yourself.**

**Note: For a correct measurement, always stand upright on the scale and make sure the electrodes are in contact with your feet.**

**Note: In order for the scale to calculate all 15 measurement indices, you must always weigh yourself barefoot. Not all body parameters can be measured if you are wearing socks. Dirty feet can also affect the results.**

**Note: Never get off the scale until the measurement is completed. When weighing yourself, your body weight is calculated first. When the measurement is completed, the numbers on the display flash twice. Then a 0 moves along the display, which completes the measurement of body proportions. When done, the results will be displayed.**

**Note: If the scale does not work properly, reset it by holding down the button that switches units from kilograms to pounds.**

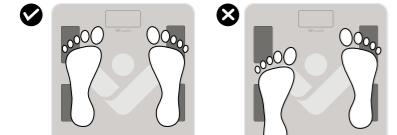Graphic representation of feet on the scale surface

### DESCRIPTION OF THE DATA ON THE DISPLAY:

**68.00**  
kg      0000 0000

The value flashes twice – weighing is complete      The number 0 moves on the display – completing the measurement of body proportions

– After weighing yourself, the scale will measure fifteen body parameters using BIA technology. Three basic parameters such as body weight, water and fat are shown directly on the display. The results are displayed 3 times in a row. You can see the other results in the application.

**MULTIPLE USERS:**  
– The scale can remember up to twelve users. To measure and display all fifteen indices, each additional user must join the current account and fill in their personal data and parameters, or register as a new user.

– After adding an account or after a new registration, the user must immediately weigh themselves so that the scale remembers them and pairs with their account.

– During weighing, CODE 2-11 appears on the top left of the screen. This CODE appears on the screen each time the user steps on.

**Note: In order for the scale to be able to recognize which user is currently weighing themselves, send the results correctly to the application and assign them to the correct account, the difference between the users must be at least 4-5 kg.**

**Note: When weighing without using the application, the results are synchronized with the application each time the scale is switched on and re-paired with the application.**

**BASIC SPECIFICATIONS:**  
Dimensions: 31 cm x 31 cm x 2.5 cm  
Weighing range: 0.2 – 180 kg  
Weighing deviation: 0.1 kg

**MANUFACTURER'S RESPONSIBILITY:**  
TrueLife products are under warranty for manufacturing defects. Eleme s.r.o. are, in the event of such defects, under obligation to repair or replace the defective part or product with a comparable part or product at its own discretion. Other alterations, conversions or further unauthorized manipulation of the product will result in the termination of the warranty. This shall also occur should continued use of defective parts, negligence or other errors in the product usage result in loss or damage. The warranty does not apply to damage resulting from standard usage, accidental or normal wear and tear. Eleme s.r.o. assumes no responsibility for any accident, injury, death, loss or other claims relating to or arising out of the use of this product. Eleme s.r.o. shall under no circumstance be made liable for any incidental or consequential damages resulting from or during the use of this product or any of its components. Repayment or replacement of spare parts and/or products may be subject to charges for transportation, handling, replacement and/or replenishment.

### SAFETY PRECAUTIONS WHILST USING THE EQUIPMENT:

1. Carefully read the enclosed safety instructions.
2. For the correct handling of the product, read the latest version of the manual available at [www.truelife.eu](http://www.truelife.eu).
3. Use only the supplied rechargeable batteries for the parent unit.
4. When using batteries in the baby unit, use the supplied non-rechargeable batteries or choose Alkaline 1.5V AAA.

Typing and printing errors are reserved.

### MANUFACTURER:

elemē s.r.o.,  
Bráškova 15, 16100 Prague 6, Czech Republic

### Explanation of terms:

Measured indices:

- BMI (Body Mass Index) is an index that is used to classify if you are underweight, overweight or various degrees of obesity. It is a calculated index from height and weight.
- AMR (Active Metabolic Rate) is an index that determines the amount of energy that the body consumes in an active state in one day.
- BMR (Basal Metabolic Rate) is an index that determines the amount of energy that the body needs at complete rest to maintain its basic functions.

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