# SonicBrush GL UV

---

**EN - USER MANUAL**

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**HR - UPUTE ZA UPOTREBU**

---

## EN - USER MANUAL


---

### BOX CONTENTS

- Body of toothbrush
- 4x Standard brush heads
- Glass
- Charging station
- Travel case
- TrueLife SonicBrush UV Sterilizer
- Manual and safety instructions

### PRODUCT DESCRIPTION


Sonic toothbrush with charging glass and travel case

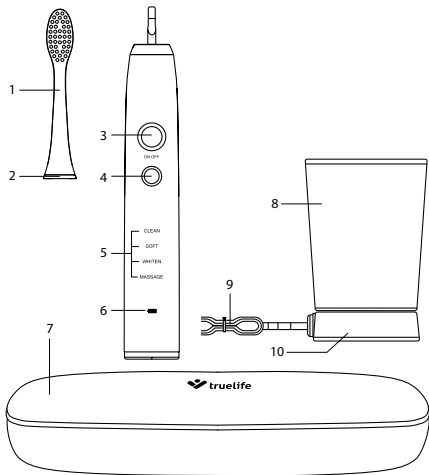- |  |                       |
|--|-----------------------|
| 1. Replaceable brush head for sonic toothbrush | 6. Charging indicator |
| 2. Separating ring                             | 7. Travel case        |
| 3. On/off button                               | 8. Glass              |
| 4. Switch mode button                          | 9. Power cord         |
| 5. Mode indicator                              | 10. Charging station  |

## UV Steriliser

1. Power cord input
2. Case cover
3. UV lamp
4. Button to unlock case cover
5. Holder for toothbrush/brush head
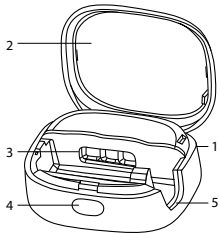
## INSTRUCTIONS FOR USE

### Sonic toothbrush

1. Place the brush head onto the body of the toothbrush.
2. Select the cleaning mode using the button. If you intend to use the same mode as in the previous use, there is no need to reset the toothbrush.
  - **CLEAN** – Regular mode for brilliantly clean teeth
  - **SOFT** – Gentle brushing for sensitive teeth and gums (may also be suitable for users less accustomed to sonic cleaning)
  - **WHITEN** – For brushing with tooth polishing
  - **MASSAGE** – For gentle gum stimulation
3. Wet the brush head and apply toothpaste to it.
4. Place the bristles of the toothbrush on the teeth, near the gums.
5. Turn the toothbrush on using the on/off button.
6. Move the bristles back and forth over your teeth and slowly move between the individual teeth. At the same time, gently push the brush towards and between the teeth
7. For best results, divide the brushing into 4 sections. Take 30 seconds for each section. After a 30-second cycle, the brush will pause and alert you to move to the next section. After 2 minutes, the brush switches itself off.
8. Rinse the brush head under running water. Wipe the body of the brush with a damp cloth and then wipe it dry with a cloth or towel .
9. Put the body of the toothbrush back into the glass or store it in the travel case.

### UV Steriliser

1. Use the button to open the lid of the UV steriliser .
2. Place the toothbrush or brush head in the UV steriliser. The bristles should face the UV lamp .
3. Close the lid. Sterilisation will begin automatically.
4. Sterilisation lasts 3 minutes. Do not open the lid during the sterilisation process.
5. Once the UV light switches off, you can open the lid and remove the toothbrush or brush head.

**Caution:** While the bristles of the toothbrush or toothbrush head may be slightly damp when placed in the steriliser, they should never be dripping wet.

**Caution:** Use the UV sterilizer only when the toothbrush or toothbrush head is facing down.

**Caution:** When transporting a toothbrush or toothbrush head inside a UV sterilizer, the bristles must be completely dry.

## CHARGING

### Sonic toothbrush

The sonic toothbrush is equipped with a contactless charging system. As a result, a direct connection to the power cord is not necessary. This makes charging safer and more convenient. To reach a full charge, the brush should be left to charge for approximately 24 hours.

1. If the charging indicator starts blinking, the toothbrush needs to be charged.
2. Connect the charging station to the adapter and then to the power supply.
3. Place the glass on the charging station.
4. Put the toothbrush (switched off) into the glass. The charging indicator lights up red while charging. When the toothbrush is fully charged, the charging indicator changes to green.

**Recommendation:** In order to ensure optimal battery life, we recommend charging the toothbrush at regular intervals (e.g., every 2–3 weeks).

**Note:** The glass for the charging station can also be used as a toothbrush stand. However, do not connect the charging base to a power source unless you are charging the toothbrush.

**Caution:** The charging station and toothbrush may become slightly warm during charging.

### UV Steriliser

If the UV steriliser needs to be charged, the blue light inside will flash 10 times. You can charge the UV steriliser using the power cord included. Follow the instructions below.

1. Connect the UV steriliser to the power cord.
2. Plug the other end of the AC adapter into the adapter and connect it to a power source. The blue indicator light will flash while charging.
3. The UV steriliser will be fully charged after approximately 3 hours. The blue light will switch off automatically once charging is complete.

**Note:** If the UV sterilizer's charge is completely empty, the LED may not light up until it is at least minimally charged.

**Note:** The adapter is not included in the package.

## MAINTENANCE

Before cleaning, make sure that the appliance is not plugged in.

### Glass

Rinse the glass regularly under running water. The glass can be cleaned in the dishwasher.

### Charging station

Wipe any dirt off the charging station using a soft cloth. Do not immerse the charging station in water.

### Toothbrush

1. Remove the brush head.
2. Wipe any dirt off the surfaces with a soft cloth and return the toothbrush to the glass. Do not immerse the whole body of the brush in water.
3. Rinse the brush head and its bristles after each use.
4. Do not wash the toothbrush or any of its components in the dishwasher.

### Travel case

Wipe the surface and inside of the travel case with a damp cloth.

### UV Steriliser

Wipe the surface and the UV steriliser holder with a damp cloth.

**Recommendation:** For hygiene reasons, we recommend changing the brush head at least every 3–4 months.

## SPECIFICATIONS

### Sonic toothbrush

Voltage input	3.7 V
Power consumption	1.7 W
Water resistance	IPX7
Battery	Li-Ion
Ambient temperature for operation	5–50 °C
Weight	538 g
Dimensions	29×300×29 mm

### UV Steriliser

Voltage input	5 V
Power consumption	0.2 W
Material	ABS
UV lamp expected life	≥5000 hours
Ambient temperature for operation	4–50 °C

UV radiation	$\geq 1000\mu\text{W}/\text{cm}^2$
Weight	33,5 g
Dimensions	60x52x25 mm

## PRODUCT LIABILITY

TrueLife products are guaranteed to be free from manufacturing defects. In the event of such defects, elem6 s.r.o. will repair or replace the defective part or product with a comparable part or product at its sole discretion. The warranty will be deemed void if the device is altered or tampered with. The warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss, or other claims related to or resulting from the use of this product. elem6 s.r.o. shall in no event be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Return or exchange of spare parts and / or products may be subject to shipping, handling, and replacement or replenishment fees.

**Printing errors and changes of information reserved.**

## SAFETY GUIDELINES WHEN OPERATING THE DEVICE

1. Read the leaflet included in the packaging carefully .
2. The latest version of the manual, including instructions for operating the device, can be downloaded at: [www.truelife.eu](http://www.truelife.eu).
3. Printing errors and changes of information reserved.

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