SonicBrush UV

Video - How to use SonicBrush UV

Video - Equipment of SonicBrush UV

EN - INSTRUCTIONS FOR USE

- Base with UV lamp
- Brush handle
- 1x Sensitive head, 1x Standard head, 1x Whiten head
- Travel case

Product Description

Brush Handle

- 1 Power button
- 2 Cleaning mode button
- 3 Mode indicator
- 4 Charging indicator

Base

- 5 UV lamp cover
- 6 Base
- 7 UV lamp
- 8 Sanitizer cover
- 9 Charging cradle
- 10 Power cord
- 11 Plug

Brush Heads

- 12 Identification ring


Charging

Before charging, make sure the brush handle is turned off. Insert the plug into the socket and place the handle into the charging cradle. When charging, the charging indicator light lights red on the handle. When fully charged the indicator turns green.

- 1 A full charge lasts for about 24 hours. Full capacity is sufficient for 20 pre-set cycles or approximately 10 days when cleaning teeth twice a day.
- 2 The brush must be charged for 24 hours:


 - Before first use
 - If it hasn't been used for 3 months
 - If the battery is completely empty

- 3 To maintain full capacity, keep the brush in the stand when idle. The stand must be plugged into a wall socket.
- 4 When charging, the brush may become hot.

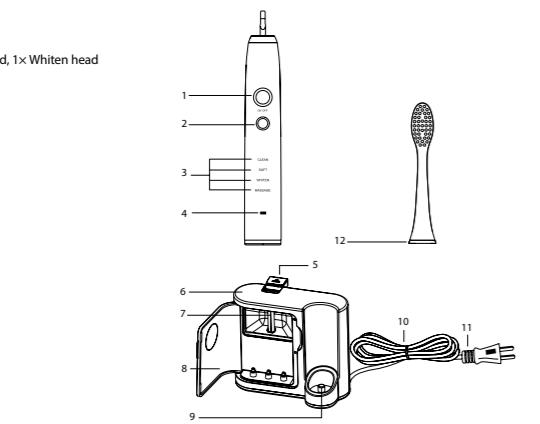Tooth cleaning

- 1 For maximum efficiency, divide the cleaning into 4 sections and give each section 30 seconds.
- 2 After 30 seconds, the toothbrush will pause briefly to alert you to move to the next section. After 2 minutes, the brush will switch itself off. Clean at least twice a day, morning and evening.
- 3 Place the brush head firmly on the metal shaft with the bristles facing the front of the handle.
- 4 Wet the head, apply toothpaste and spread it over the teeth.
- 5 Turn the toothbrush on with the ON/OFF button. To select another mode, press the button below.
- 6 In each section, start by cleaning the outside of the teeth, then the inside of the teeth, and finally clean the surfaces. Hold the toothbrush tilted so that the bristles gently touch the teeth and gum margins at the same time. When brushing, do not push the brush and do not make circular movements as with a manual toothbrush, but move slowly from one tooth to the next.
- 7 Rinse the brush head and body under running water and then wipe with a cloth or towel and place back in the base.

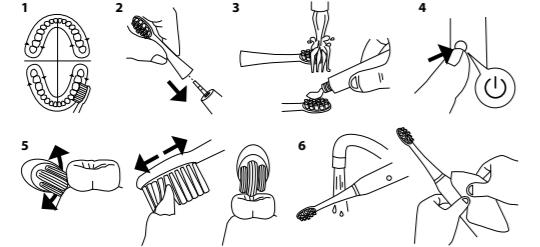Cleaning modes

CLEAN - For routine cleaning and plaque removal

SOFT - For cleaning sensitive teeth and gums

WHITEN - For removing stubborn plaque and whitening teeth

MASAGE - For stimulating the gums and improving their health

UV Lamp

The base includes a UV lamp that disinfects the cleaning head after each use. The UV lamp can remove 99.9 % of common bacteria such as E. Coli, Salmonella, Candida Albicans and the H. Simplex or H1N1 viruses.

- Only disinfect the original head with the UV lamp.

- The cleaning cycle lasts around 10 minutes and when finished the lamp switches itself off.

- The cycle starts automatically when the sanitizer cover is closed.

1. Clean and dry the head.

2. Open the UV lamp cover.

3. Place the head in the rack with the bristles pointing towards the UV lamp.

4. Make sure the charging base is plugged in.

5. Close the sanitizer cover and the disinfection process automatically begins.

- The blue illumination of the UV lamp indicates that the lamp is on and the cleaning cycle is in progress. The light might not be noticeable in brightly lit spaces.

- After the cleaning cycle is finished, the lamp automatically turns off and the light goes out.

Maintenance

Before cleaning, make sure the device is not plugged in.

Base

Wipe the dirt from the charging cradle and UV lamp with a soft cloth.

Brush

1. Remove the brush head.

2. Wipe the dirt off all surfaces with a soft cloth and put the brush back on the charging base. Do not immerse the handle completely in water.

3. Rinse the brush head and bristles after each use.

4. Do no wash the brush or any of its parts in the dishwasher.

Recommendation: Change the heads after 3 months or when the blue bristles turn white.

Troubleshooting

Problem	Possible Cause	Solution
Brush does not work	Brush has never been turned on or has not been used for more than 3 months	Allow brush to charge for at least 24 hours

Weak or intermittent vibration

Brush is pressed too hard on the teeth

Reduce the amount of pressure on the teeth

Brush only works for a short amount of time or switches off immediately when disconnected from the charging base

Battery is at the end of its life cycle

Contact customer support

Sanitizer cover is not closed

Close cover

Base not plugged in

Plug the base into the wall socket

UV lamp is damaged

Contact customer support

- If the UV lamp shatters, do not touch it and leave it for 30 minutes after which time it is safe to clear the shards and contents.

Responsibility of the Manufacturer

TrueLife products are subject to a warranty that covers manufacturing defects. In the event of such defects, the company elem6 s.r.o. is obliged under warranty to either repair or replace the complete product as we may choose. If any conversion, alteration or repair is undertaken by an unauthorized person, the guarantee becomes void. The warranty does not cover any damage resulting from improper use, normal wear or caused accidentally. The company elem6 s.r.o. assumes no responsibility for any accidents, injuries, death, loss or other claims relating to or arising from the use of this product. The company elem6 s.r.o. in no event shall be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Repayment or replacement of parts and/or products may be subject to shipping, handling, replacement and/or shipping charges.

Printing error and information changes reserved.

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